

Pure brilliance

New Keeler *Jazz* Diagnostic Set



Keeler

# Keeler


Keeler have been manufacturing Ophthalmoscopes and Otoscopes for almost 100 years.

We only supply high quality precision optics providing the best possible views.

Our latest introduction to our hand held range is no exception. Pure Brilliance is what you get from our new Jazz range of instruments...read on to discover the new Keeler Jazz range.

## Contents

Introduction .....	p3
Jazz LED Ophthalmoscope .....	p4
Jazz Fibre Optic LED Otoscope .....	p5
LED Illumination and Batteries .....	p6
Part Numbers and Accessories .....	p7
Contact details .....	p8

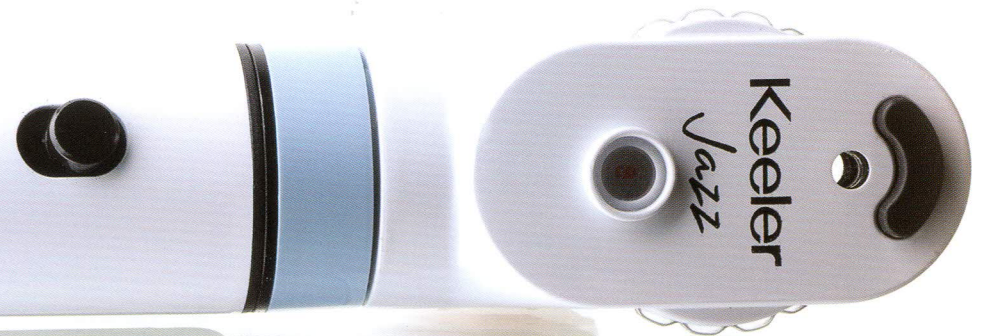






# Pure brilliance

That is exactly what you achieve every time you use your Keeler Jazz instruments.

Designed to fit neatly into your pocket, the new Jazz still encompasses all of the features you would expect from a Keeler instrument.

And why not change the look of your Jazz instrument by simply swapping over the colour rings...choose from pink, blue or black.

Supplied with 3 interchangeable colour rings

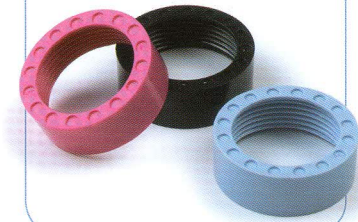




# New Jazz LED Ophthalmoscope

## High Performance Optics


Our renowned reputation for quality optics ensures you get the best possible views without compromise.

## Optical Pathway


Corneal and iris reflexes are avoided by separation of the parallel observation and illumination beams (Gullstrand's principle). The illuminated retinal area is fully visible – even with narrow pupils. Optimum conditions are achieved for ophthalmological examinations.

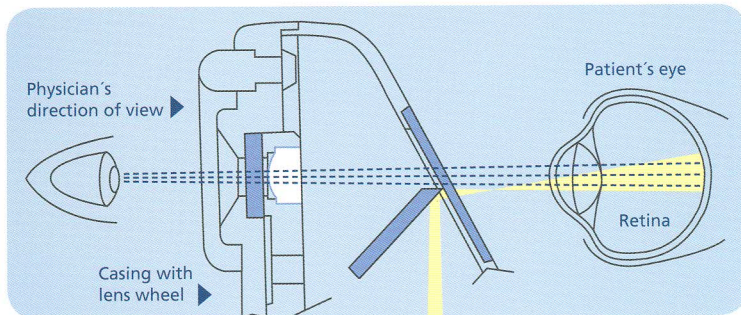

Healthy fundus

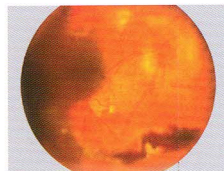

Pathological fundus

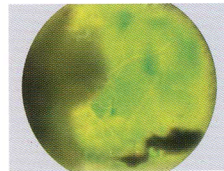

Fundus with red free filter

## Apertures and Filters



- Brilliant LED illumination provides crisp, clear views for easier diagnosis
- Trusted Keeler optics providing you with the clearest, crispest views of the retina
- Lens powers +/- 1, 2, 3, 4, 6, 8, 10, 15 and 20 Dioptres
- Smooth aperture wheel for precise selection of apertures
- Protects your spectacles



lighter



brighter



greener



long lasting


- 1 On/Off switch
- 2 Dioptre display
- 3 Dioptre adjustment wheel
- 4 Aperture wheel
- 5 Interchangeable colour rings
- 6 Spectacle protection
- 7 Handle/battery compartment





## New Jazz Fibre Optic LED Otoscope

- Pure white uniform illumination shows colours as they really are for more efficient diagnosis plus longer life battery and no more replacement of bulbs
- Sealing precision optics with a triple magnification for optimum views every time
- Can be used for performance of pneumatic tests (need connector and insufflator)
- Specula available in 5 sizes, 2mm, 2.5mm, 3mm, 4mm and 5mm

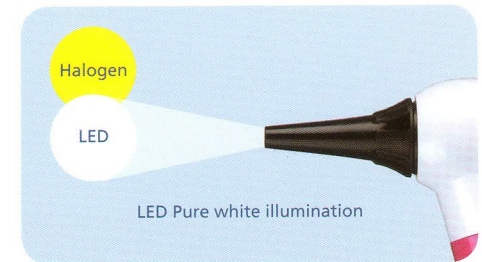

lighter



brighter



greener



long lasting

- 1 On/Off switch
- 2 3x magnifying pivoting lens
- 3 Pneumatic test connection
- 4 Fibre optics
- 5 Interchangeable colour rings
- 6 Ear speculum
- 7 Handle/battery compartment






# LED Illumination and Batteries

## Green efficiency

The new Keeler Jazz range of Ophthalmoscopes and Otoscopes with LED technology provide you with bright white illumination, as well as longer life batteries and longer illumination life.

You can achieve up to 20,000 hours illumination and up to 50 hours battery life with two conventional AA batteries.

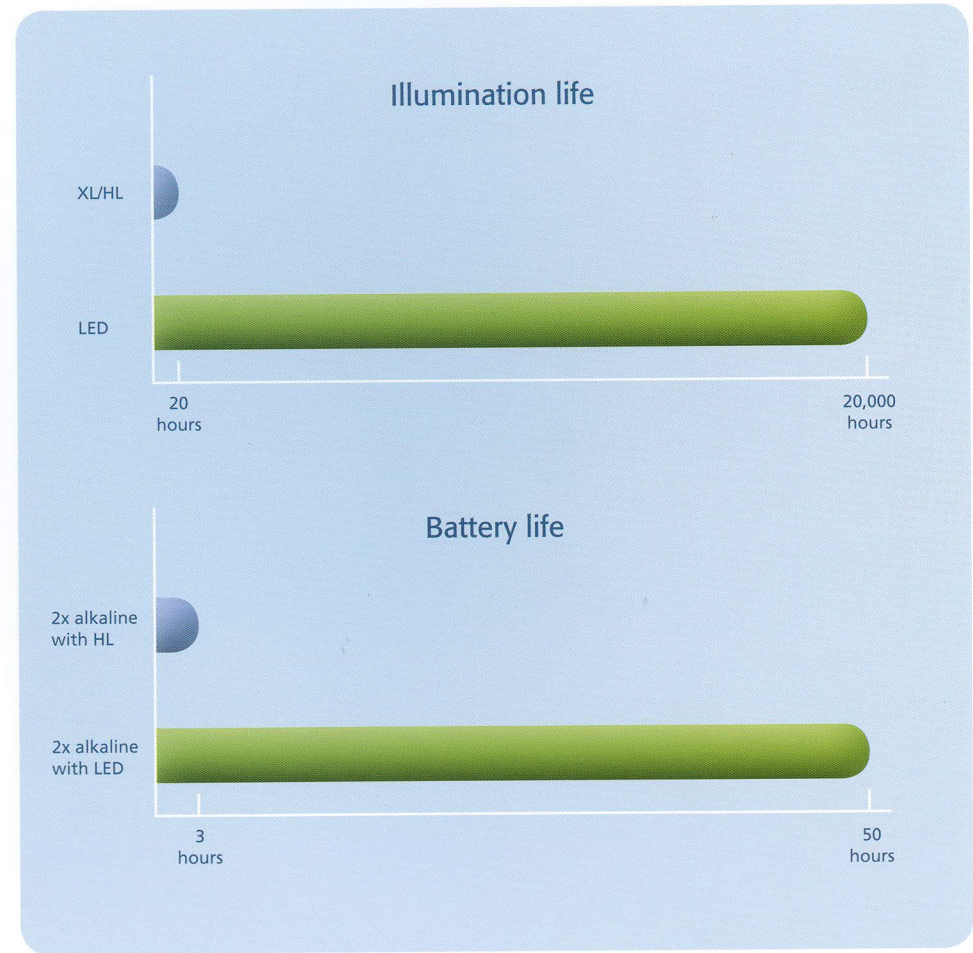



## Part Numbers and Accessories



**Jazz Ophthalmoscope  
in Single Pouch**

Part Number **1102-P-1092**


**Jazz LED  
Ophthalmoscope  
in Single Pouch**

Part Number **1102-P-1105**



**Jazz LED Otoscope in  
Single Pouch**

Part Number **1514-P-1160**

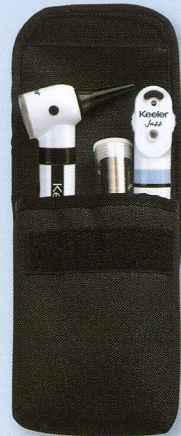

**Jazz Diagnostic Set  
in Double Pouch**

Part Number **1702-P-1088**

**Jazz LED Diagnostic  
Set in Double Pouch**

Part Number **1702-P-1096**




**Ophthalmoscope Bulb**

Part Number **1011-P-5530**



**Ophthalmoscope LED**

Part Number **1011-P-5522**

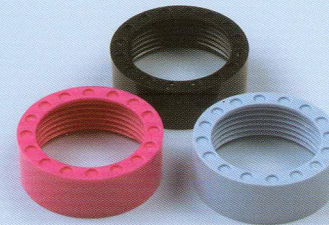

**Spare Coloured Rings**

Part Numbers

**EP39-37051** Black Ring

**EP39-37342** Pink Ring

**EP39-37350** Blue/Grey Ring



**Disposable Specula**

Part Numbers

**1514-P-7056** 2mm Pack of 10

**1514-P-7044** 2.5mm Pack of 10

**1514-P-7052** 3mm Pack of 10

**1514-P-7060** 4mm Pack of 10

**1514-P-7079** 5mm Pack of 10

**1514-P-7087** 2mm Pack of 100

**1514-P-7095** 2.5mm Pack of 100

**1514-P-7108** 3mm Pack of 100

**1514-P-7116** 4mm Pack of 100

**1514-P-7124** 5mm Pack of 100