



Vitreo

KUHN 45° POWER FORCEPS



Continuing its leading role in the development of ophthalmic instrumentation, Vitreo has launched a new, unique, all-purpose forceps for standard and heavy duty membrane peeling. This new Kuhn 45° Power Forceps has been developed in collaboration with Professor Ferenc Kuhn.


- The Kuhn 45° Power Forceps is designed for peeling proliferative intravitreal, epi and sub-retinal membranes, lens capsule, iris, small IOFB, fibrinous membranes and more. The forceps is also suitable for peeling standard ILM and ERM membranes and attached posterior hyaloids.
- Its unique tip design enables exceptional visualization while grasping the target tissue, with a 45° crocodile design inner surface that enables selective and strong grasping of the target tissue, regardless of its rigidity, size, orientation, or surface texture.



“I have not come across any task in the eye that I am not able to do with the Kuhn 45° Power Forceps.”

Professor Ferenc Kuhn

- The VR surgeon can easily check internal alignment of the tip of the forceps during surgery via special marking on the instrument close to the tip. The forceps also features an ergonomically designed handle for improved intraocular control and precision during surgical maneuvers.

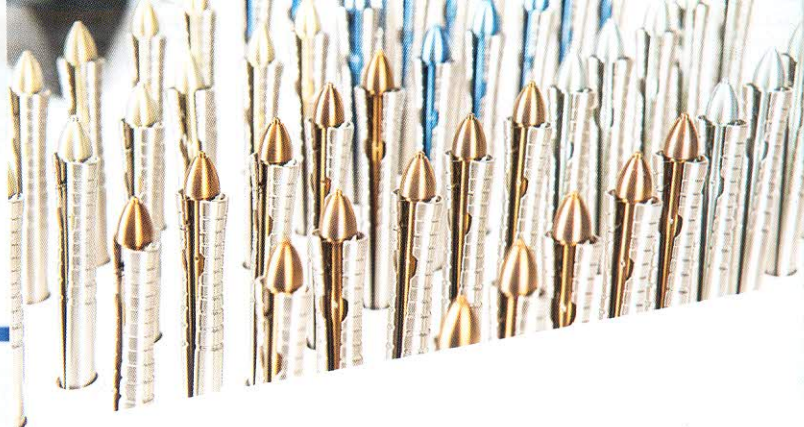

Available in 23G and 25G in both reusable and disposable versions. The reusable forceps will be delivered including the unique cleaning and sterilization holder.



The Vitreq cleaning and sterilization holder.

23 Gauge



Packaged box/5, Sterile

25 Gauge



Packaged box/5, Sterile

About Vitreq

Vitrex creates instruments for the ophthalmic society with a special focus on vitreoretinal surgery. Our products are perfectly matched to your needs. Fluid design that is achieved through working hand-in-hand with leading, innovative surgeons from all over the world: That is our vision.

With many years experience in design, production and marketing of high quality products, Vitreq develops reliable and safe, next generation solutions that enhance performance and support advanced surgical techniques.



Au **DISTRIBUCI ZAJIŠTUJE:**



SPRIT MEDICAL spol. s r.o.
Brno, Sadařská 495/1, PŠČ: 624 00
Tel.: 548530828, fax: 548538897
E-mail: objednavky@spiritmedical.cz
[http:// www.spiritmedical.cz](http://www.spiritmedical.cz)

**For international
please contact
Vitrex Headquarters**
Seggelant-Noord 2
3237 MG Vierpolders
The Netherlands
T: +31 181 745090
F: +31 181 478583
E: info@vitreq.com

**In USA please contact
Vitrex USA Inc.**
2 Marshall Rd. Unit G
Kingston, NH 03848
United States
T: +1 603 347 5590
F: +1 603 347 5583
E: info@vitrequsa.com

**In Germany please
contact
Vitrex Deutschland**
Prinzenallee 7
D-40549 Düsseldorf
Deutschland
T: +49 7251 36 97 123
F: +49 7251 36 97 124
E: info@vitreq.de

www.vitreq.com

